


Today's technology guarantees high quality chips

Safe and easy way to improved bone chips

The Noviomagus Femoral Head Reamer Set is designed to fixate different types of bone material for either cartilage removal, shaping or splitting. The set comprises a base plate, a detachable clamp, a screw handle and a splash shield. Different sizes of inverted reamers are available to remove cartilage from a femoral head. Removing cartilage improves remodelling of the bone graft and enhances osseointegration of the graft.

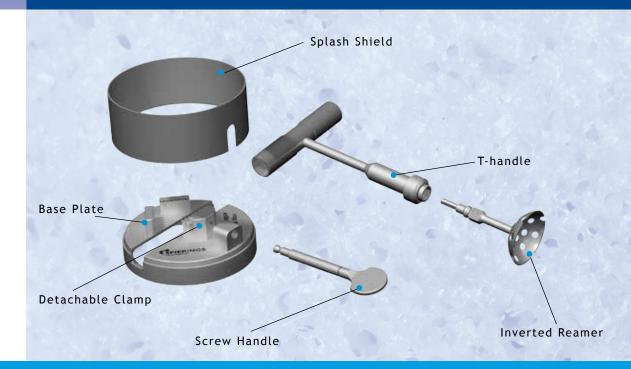

Safe and solid fixation of bone material

Removing cartilage from the bone requires a safe and solid fixation of the bone material. The bone vice is designed for this specific purpose. An adjustable clamp guarantees accomodation of every size and shape of bone material. The sharp teeth make sure the bone material is well fixed. A slot in the bottom of the bone vice makes it easy to split the bone by means of a saw or osteotome. The bone vice is equipped with a reusable splash shield, designed to prevent the spread of the removed tissue and contamination of the reaming area.





Reams fresh frozen bone

The inverted reamers are made of a special hardened stainless steel alloy, resulting in a high quality and a long life-span product. When fitted on a surgical power tool, the inverted reamers can remove the cartilage from a femoral head within 30 seconds. Thawed and fresh frozen femoral heads are reamed with equal ease. This avoids the risk of defrosting too much bone prior to surgery. The shape of the reamers make sure that only cartilage will be removed. This results in time- and money-saving efficiency during surgery.




Within 30 seconds!

Improve the quality of your revision

For the quality and longevity of graft reconstructions it is necessary to remove the cartilage. Cartilage does not show any signs of incorporation within the grafted area. Research shows that cartilage affects the stability of the graft and increases the risk of loosening in the long-term¹. Using the inverted reamers to remove cartilage is a success-guaranteed

method. The Noviomagus Inverted Reamers are available in 10 diameters (stepwise increasing from Ø40 mm to Ø58 mm, with increments of 2 mm) covering all sizes of femoral heads. They are available with quick-release couplings that fit the standard A0-coupling, Zimmer-Stryker-Hall power tools and the specially designed manual T-handle.

 S. van der Donk, P. Buma, T.J.J.H. Slooff, J.W.M. Gardeniers, B.W. Schreurs: Incorporation of Morselized Bone Grafts: a study of 24 acetabular biopsy specimens. Clinical Orthopaedics and Related Research (2002) 396: 131-141.

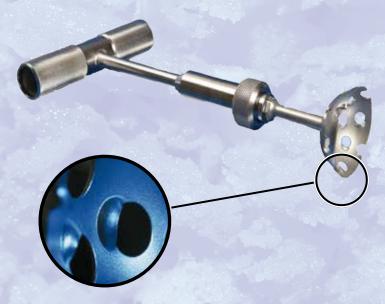

Inverted Reamers only remove cartilage, no precious bone material

Advanced benefits

- Saving time in surgery
- User-friendly design
- Focus on safe cartilage removal
- Improving bone graft quality
- Fast cartilage removal of fresh frozen bone
- Inverted Reamer diameter ranging from ø40 to ø58 mm
- 3 Years warranty for regular hospital use
- Reamers suitable for both manual and power tool handing


Ordering information

Noviomagus

Femoral Head Reamer Set

Noviomagus Bone Vice:

Noviomagus Bone Vice:	
Base Plate	SBM-1109-14
Detachable Clamp	SBM-0503-15
Screw Handle	SBM-0503-16
Splash Shield	SBM-1109-17
Sterilization Tray (NFHRS)	SBM-1109-45
T-handle AO-coupling	SBM-1109-43
Inverted Reamer-AO Ø 40mm	SBM-1303-62
Inverted Reamer-AO Ø 42mm	SBM-0503-18
Inverted Reamer-AO Ø 44mm	SBM-0904-21
Inverted Reamer-AO Ø 46mm	SBM-0904-22
Inverted Reamer-AO Ø 48mm	SBM-0503-19
Inverted Reamer-AO Ø 50mm	SBM-0904-23
Inverted Reamer-AO Ø 52mm	SBM-0904-24
Inverted Reamer-AO Ø 54mm	SBM-0503-20
Inverted Reamer-AO Ø 56mm	SBM-0904-25
Inverted Reamer-AO Ø 58mm	SBM-0904-26
Inverted Reamer Z-S Hall Ø 42mm	SBM-1207-50
Inverted Reamer Z-S Hall Ø 44mm	SBM-1207-51
Inverted Reamer Z-S Hall Ø 46mm	SBM-1207-52
Inverted Reamer Z-S Hall Ø 48mm	SBM-1207-53
Inverted Reamer Z-S Hall Ø 50mm	SBM-1207-54
Inverted Reamer Z-S Hall Ø 52mm	SBM-1207-55
Inverted Reamer Z-S Hall Ø 54mm	SBM-1207-56
Inverted Reamer Z-S Hall Ø 56mm	SBM-1207-57
Inverted Reamer Z-S Hall Ø 58mm	SBM-1207-58

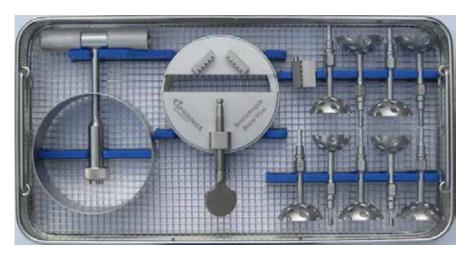

The Noviomagus Femoral Head Reamer Set is CE marked and FDA approved. Other Noviomagus products manufactured by Spierings
Orthopaedics B.V. are: Noviomagus Impactor Set,
Noviomagus Mini Mill, Noviomagus Bone Mill and the Noviomagus Revision Meshes.

No rights can be derived from the information in this brochure.

Copyright: Spierings Orthopaedics B.V.











Manufactured by: Spierings Orthopaedics B.V. Madoerastraat 24 6524 LH Nijmegen The Netherlands

T: +31 -(0)24-3501603 F: +31 -(0)24-3501604 E: info@spierings.biz

I: www.spierings.biz